
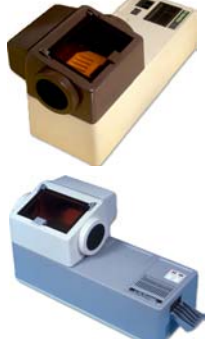



Peri-Pro I, II & III Model Comparison Sheet

Model	Peri-Pro Original	Peri-Pro II	Peri-Pro III
Picture	 <p>1977-1989</p> <p>1989-present</p>	 <p>1982-1989</p> <p>1989-1997</p>	 <p>Only 1 model color</p>
Manufacture Date	1977 - Present Day	1982-1997	1997 - Present day
Processing Time	7 minutes 30 seconds	6 minutes 30 seconds	5 minutes
Chemistry Temp	Room Temp 70 - 82 ° F	74 - 76 ° F	75 ° F
Chemistry Heater	NO	Heater Bar in both the chemistry tanks	Heater Plate underneath the chemistry tanks
Temperature Sensor	NO	Thermistors located inside the Heater Bar	Thermsitor probe in the developer tank
Film Dryer	Dryer Heater Element 350 watt 45 ohms	Dryer Heater Element 350 watt 45 ohms	Dryer Heater Element 370 watt 35 ohms
Thermal Fuse for Film Dryer	112 ° C / 233.6 ° F trip	112 ° C / 233.6 ° F trip	98 ° C / 208.4 ° F trip
Main Power Fuse	5 Amp	5 Amp	5 Amp
Power Switch	Push Button Serial # 14999 & ↓ Rocker Style Serial # 15000 & ↑	White Switch S/N 15142 & ↓ Rocker w/ light S/N 15143 & ↑	Rocker Style Switch Amber Light
Process Switch	No Process Switch	White Switch S/N 15142 & ↓ Rocker w/ light S/N 15143 & ↑	Rocker Style Switch / Green ready light
Main Drive Gear	30 tooth gear #304	36 tooth gear #305	66 tooth gear
Transport Timing	Timing Gear dot in 10:00 o'clock position on clock face	Timing Gear dot in 10:00 o'clock position on clock face	Timing Gear dot in 1:30 o'clock position on clock face
Maintenance Door	No, cover has to be removed	No, cover has to be removed	Yes
Water Drain Valve	No, cover has to be removed	No, cover has to be removed	Yes